

Medical Garment and Personal Protective Equipment





About Us

Medtecs International Corporation Limited ("Medtecs") is a leading supplier and distributor of personal protective equipment ("PPE") and provider of hospital logistics services with over 30 years of experience in the manufacturing of PPE and workwear.

Medtecs and its subsidiaries (collectively, the "Group") commenced operations in 1989 and have since established a strong presence in the United States, Europe and Asia Pacific. The Group has offices and facilities spanning across Asia including Singapore, Taiwan, the Philippines, the People's Republic of China, Cambodia, and the United States of America for expansion into the global domestic PPE and medical consumables market. Currently, the Group's headquarters is located in Taipei, Taiwan.

Medtecs is dual-listed on Singapore Exchange Limited (SGX: 546) and Taiwan Stock Exchange (TWSE: 9103).

The Group is dedicated to safeguarding the health and safety of people worldwide.



Medtecs (Cambodia) Corporation Limited

Disposable PPE, Medical Garment, Workwear, Apparel



RMKH Glove (Cambodia) Co., Ltd.

Nitrile Glove



Medtecs (MSEZ) Corporation


Disposable PPE



Medtecs International Corporation Limited (Philippine branch)

Face Mask, Medical Garment, Uniform

Our Strengths

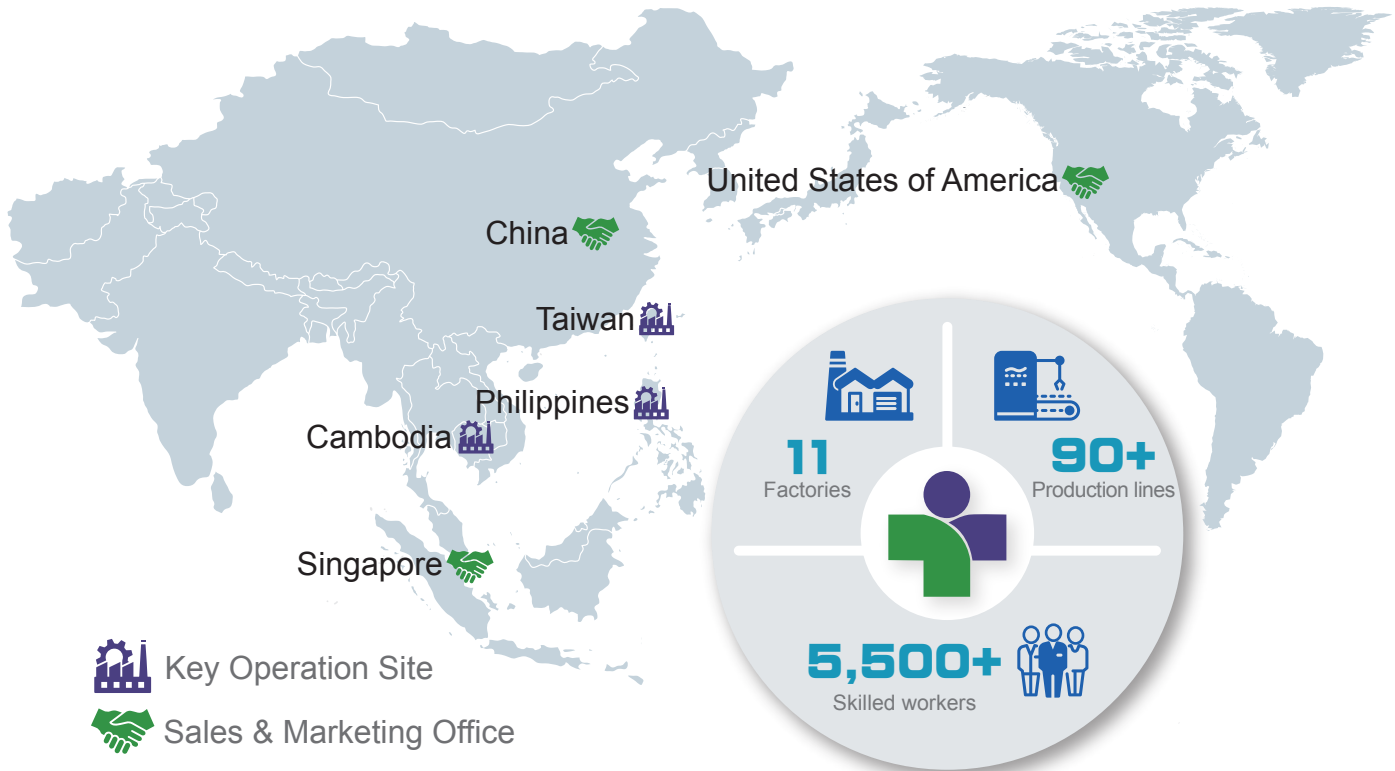

World's #1 Capacity in Seam Taping Processing

We have the world's largest capacity in seam taping processing, meeting the growing demand for personal protective equipment.

Designated Supplier for Leading Brands

We have built long-standing relationships with our clients through decades of experience in manufacturing and supplying high-quality personal protective equipment.

Agile Capacity and On-Time Delivery

Highly scalable and flexible manufacturing capacity with strong distribution channels to fulfill diverse needs and guarantee on-time delivery.

Non-China Based Manufacturing Sites

We maintain diversified manufacturing facilities and production lines spanning across Asia, including the Philippines, Cambodia, and Taiwan. In particular, Cambodia and Taiwan are designated countries under US Trade Agreement Act (TAA).

Versatile Product Range Disposable & Reusable

With a comprehensive product line of disposable and reusable medical garments, we can meet all protection needs and customer demands.



Table of Contents

Reusable Medical Garment

For Medical Professionals

| | |
|---|---|
| Lab Coat _____ | 5 |
| Nurse Uniform _____ | 6 |
| Reusable Isolation Gown _____ | 6 |
| Reversible Scrub Suit _____ | 7 |
| Scrub Suit _____ | 7 |
| Scrub Suit (4-Way Stretch Fabric) _____ | 8 |

For Patients

| | |
|----------------------------|----|
| Infant Shirt _____ | 9 |
| Nursing Gown _____ | 9 |
| Pajamas _____ | 10 |
| Regular Patient Gown _____ | 10 |

Accessories

| | |
|-----------------------|----|
| Blanket _____ | 11 |
| Laundry Bag _____ | 11 |
| OR Sheet _____ | 12 |
| Sheet & Pillow _____ | 12 |
| Thermal Blanket _____ | 12 |
| Underpad _____ | 13 |

Disposable Personal Protective Equipment

For Medical Professionals

| | |
|--|----|
| Coverall (With or Without Tape) _____ | 15 |
| Dental Gown _____ | 16 |
| High Risk Isolation Gown with Tape _____ | 16 |
| Isolation Gown _____ | 16 |
| Lab Coat _____ | 17 |
| Scrub Suit _____ | 17 |
| Surgical Gown _____ | 18 |
| Warm Up Jacket _____ | 18 |

Accessories

| | |
|----------------------------------|----|
| Boot Cover _____ | 19 |
| Bouffant Cap _____ | 19 |
| Glove _____ | 19 |
| Medical/Surgical Face Mask _____ | 20 |
| Shoe Cover _____ | 20 |
| Stretcher Sheet _____ | 20 |
| Transfer Sheet _____ | 20 |



Reusable Medical Garment

For Medical Professionals



Lab Coat

- Lapel collar.

Customization

Style:

Snap front closure or with buttons.
With or without knit cuffs.

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free

Nurse Uniform

- Stylish, functional and comfortable.

Customization

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester
72% Polyester/21% Rayon/7% Spandex

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free



Reusable Isolation Gown

- Passes EN13759 and AAMI PB-70.
- AAMI level 2/AAMI level 3 standard.
- Knit cuff style.
- Tie closure to secure overlap and prevent gaping.
- Withstands 100 home washes or 25 industrial washings.

Customization

Material:

Polyester
Cotton
T/C
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free



For Medical Professionals

Reversible Scrub Suit

- V-Neck scrub top.
- Reversible top and pants.

Customization

Style:

Pants featuring a color-coded drawstring or elastic waistband.

Material:

- Polyester
- Cotton
- T/C
- CVC
- 100% BCI Cotton
- 100% ECO Polyester
- 72% Polyester/21% Rayon/7% Spandex

Finishing:

- Anti-bacterial
- Anti-static
- Water repellency
- Wrinkle free



Scrub Suit

Customization

Style:

V-Neck/U-Neck with buttons.
Pants featuring a color-coded drawstring or elastic waistband.

Material:

- Polyester
- Cotton
- T/C
- CVC
- 100% BCI Cotton
- 100% ECO Polyester
- 72% Polyester/21% Rayon/7% Spandex

Finishing:

- Anti-bacterial
- Anti-static
- Water repellency
- Wrinkle free



Scrub Suit (4-Way Stretch Fabric)

- 4-way stretch for perfect movement.
- Anti-static fabric.
- Wrinkle-resistant with odor control.
- Hidden interior pocket to store your valuable items.
- Locker loop included.
- Machine wash & tumble dry under 70°C.

Customization

Material:

72% Polyester/21% Rayon/7% Spandex

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free

For Patients

Infant Shirt

Customization

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Wrinkle free



Nursing Gown

- Comfortable and practical.
- 2 nursing slits.
- 4 snaps on IV sleeves.
- Heavyweight construction.

Customization

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free



Pajamas

Customization

Style:

With lapel collar or V-neck.
Long/Short sleeves.

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free

Regular Patient Gown

- Complete coverage and modesty design.
- 100% polyester twill tape adds durability.

Customization

Style:

Side-tie back
Side-overlap back
Straight back

Material:

Polyester
Cotton
T/C
CVC
100% BCI Cotton
100% ECO Polyester

Finishing:

Anti-bacterial
Anti-static
Water repellency
Wrinkle free



Accessories



Blanket

Material:

- 100% Polyester
- 55% Cotton/45% Polyester

Laundry Bag

- Buckle closure.
- Traditional laundry bag style.

Material:

- 100% Polyester





OR Sheet

Material:

100% Cotton

Sheet & Pillow

T140: 114 gsm \pm 5%

T180: 105 gsm \pm 5%

T250: 125 gsm \pm 5%

Material:

55% Cotton/45% Polyester

Thermal Blanket



- FR Test Standard: BS EN ISO12952.
- 2 sides or 4 sides hem.

Material:

100% Polyester

55% Cotton/45% Polyester



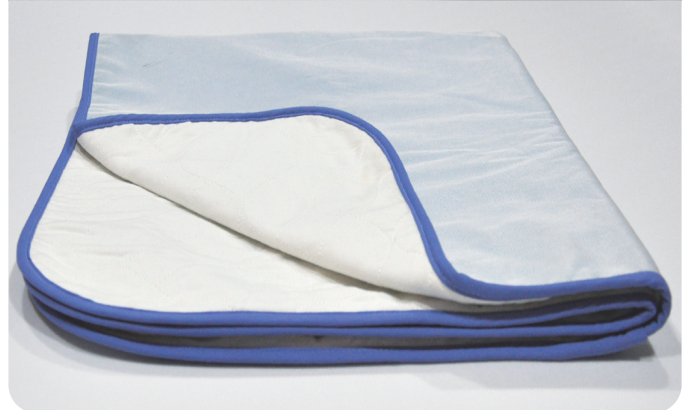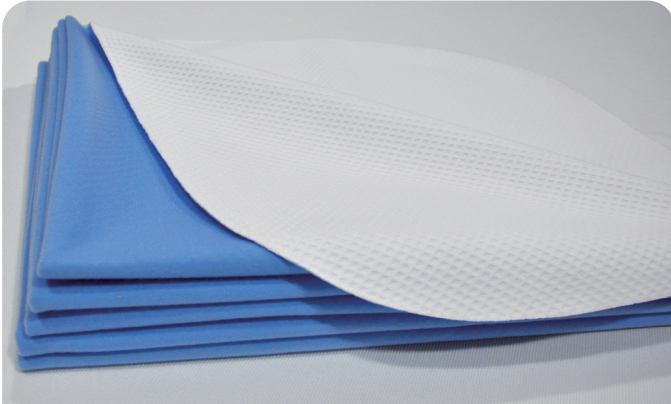

Underpad

Customization

Material:

100% Polyester

Surface:

80% Cotton/20% Polyester

Soaker:

70% Polyester/30% Rayon

Barrier:

100% PVC



Disposable Personal Protective Equipment

For Medical Professionals



Coverall (With or Without Tape)

- PPSB laminated PE.
- AAMI level 4 standard/ASTM F1670 standard/ASTM F1671 standard.
- PPE Cat. III Type 5-B/6-B.
- Breathable and impervious.
- Hood design provides overall protection.
- With elastic wrist and open ankle.

Material:

PPSB laminated with PE film

Fabric weight:

63 gsm

Dental Gown

- Made of breathable and fluid-resistant material.
- Tie at waist for wearer to adjust fit.
- Knit cuff style.

Material:
SMS/SMMS



High Risk Isolation Gown with Tape

- PPSB coated PE materials provide excellent protection from liquid-borne microbes.
- Ultra-sonic welded or sealed seams with tape design provides high performance on hydrostatic pressure and fluids impact.
- Tie at the neck and waist for secure protection.
- AAMI level 4 standard/ASTM F1670 standard/ASTM F1671 standard.

Material:
PP/PE

Fabric weight:
36-40 gsm



Isolation Gown

- Breathable SMS/SMMS providing fluid resistance.
- Tie at waist.
- AAMI level 2 standard/ASTM F1670 standard.

Material:
SMS/SMMS

Fabric weight:
25 gsm



Lab Coat

- Made of breathable material, this garment is not only fluid-repellent but also features anti-static properties.
- Designed with traditional collar and knit cuff.

Material:

PPSB or SMS/SMMS



Scrub Suit

- Made of breathable and fluid-resistant material.
- Special treatment for antistatic function.
- With one chest pocket and two pockets at waist.

Material:

PPSB or SMS/SMMS



Surgical Gown

- Made of breathable and fluid-resistant material.
- Features alcohol-repellent and anti-static treatments.
- Ultra-sonic welded and heat sealed seam.
- 4 ties at waist to prevent cross contamination.
- Knit cuff style.
- AAMI level 2/AAMI level 3 standard.

Material:

SMS/SMMS

Fabric weight:

40-50 gsm



Warm Up Jacket

- Made of soft non-woven material.
- Knit cuff style.
- Front snap closure.

Material:

PPSB or SMS/SMMS

Accessories



Boot Cover

- Made of PPSB coated with PE, non-breathable but fluid-impervious.
- Ultra-sonic welded seams to prevent blood and fluids from penetrating through pinholes.
- Anti-slip treatment on bottom provides additional traction.
- With adjustable strap.

Material:

PPSB with PE coating

Bouffant Cap

- Made of breathable PPSB.
- Puffy shape, typically secured around the head with elastics.
- With or without anti-static.

Material:

PPSB

Packaging:

100 ea/bx, 1000 ea/ctn



Glove

- Highly elastic and durable.
- Comfortable fit.
- Odorless and hypoallergenic.

Material:

Nitrile Butadiene Rubber (NBR)

Cuff:

Beaded, 9 inches ± 1 inch.

Color:

Blue/Cool Blue/Cobalt Blue/Sky Blue/Black/White

Size:

XS-XL

Packaging:

200 pcs/bx for size XS-L

180 pcs/bx for size XL

10 bxs/ctn



Medical/Surgical Face Mask

- Premium 3-ply mask.
- 99% Bacteria Filtration Efficiency.
- 99% Virus Filtration Efficiency.
- Over 80% Particulate Filtration Efficiency.
- Synthetic blood penetration resistance test passed.
- Pleated design.

Material:

Outer Layer-PPSB
 Middle Layer-Meltblown
 Inner Layer-PPSB



Shoe Cover

- Non-skid, heavy duty, fluid-resistant, low lint.

Material:

PPSB/PPSB with PE coating



Stretcher Sheet

- Fluid repellent.
- Fits standard and large stretchers/gurneys.
- Elastic ends with deep corner pockets.
- Reduces laundry/linen costs.

Material:

PPSB/SMS

Fabric weight:

50 gsm

Color:

White

Size:

Regular, XL



Transfer Sheet

- Saves laundry time and costs.
- Helps to prevent cross contamination.

Material:

PPSB

Fabric weight:

90 gsm

Color:

Navy Blue

Size:

Regular, XL

Taiwan

Medtecs (Taiwan) Corporation

11F, No.9, SongGao Rd., Xinyi Dist., Taipei City 110, Taiwan
+886-2-2739-2222

Singapore

Medtecs (Asia Pacific) Pte. Ltd.

138 Cecil Street, #13-02 Cecil Court, Singapore 069538
+65-6534-9293

Philippines

Medtecs International Corporation Limited – Philippine Branch

22F, The World Center Building, #330 Sen. Gil Puyat Avenue Bel-Air,
Makati City, Philippines
+632-8555-5621

Cambodia

Medtecs (Cambodia) Corporation Limited

No. 52, Street 606 Corner Street 311, Phum 8,
Sangkat Boeng Kak 2, Khan Tuol Kouk,
Phnom Penh, Cambodia, 120408
+855-23-866-659

China

Hangzhou Jincheng Medical Supplies Manufacture Co., Ltd.

202 Zhangshan Road, Yuhang District, Hangzhou, China 311107
+86-571-8639-6888

Bermuda

Medtecs International Corporation Limited

Clarendon House 2 Church Street, Hamilton HM11, Bermuda
+63-2-8817-9000

United States

Medtecs USA Corporation

1390 Market Street, Suite 200, San Francisco, California 94102, USA



Medtecs International Corporation Limited

<https://business.medtecs.com>

Email: sales@medtecs.com

